


LAB 16. RAY OPTICS: HOLLOW LENS

AP PHYSICS II

Driving Question | Objective

Is a Convex Lens always convergent? Is a Concave Lens always divergent? What are some factors which would affect the behavior of a lens?

You will attempt to build an understanding on the connection between indices of refraction and geometric configuration of a lens in its behavior to converge or diverge light.

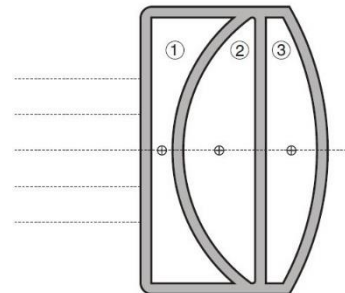

Design and Conduct Your Experiment


It is your group's responsibility to design and conduct an experiment whose data will support your answer to the driving question above. Use the answers to the pre-lab questions below to help guide your experiment design. After you have answered the pre-lab questions, write an outline of the equipment setup and procedure you will use to collect data, identifying the steps in sequence and the points at which each piece of equipment will be used.

Materials and Equipment

- Light Source (Ray Side)
- Hollow Lens
- Lens Box
- Water
- Pipette
- 150 ml Beaker

Experimental Design

1. Set the light source to output 5 parallel rays.
2. Place the light source on the left side of the box and shine the incident rays into the box. Assure that the incident rays are perpendicular to the flat side of the hollow lens.

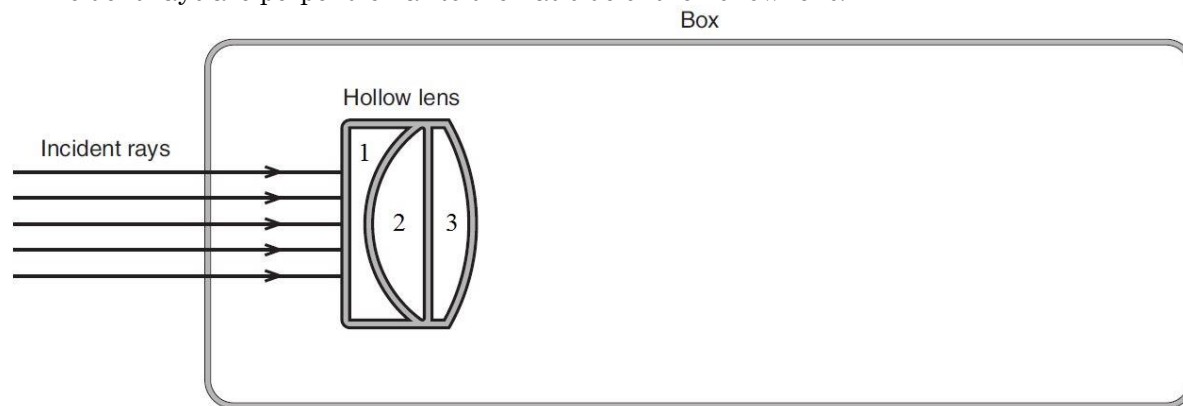

3. Is there significant refraction occurring? If so, is it convergent or divergent?

There is no significant refraction

4. We will be adjusting the indices of refraction of each portion of the hollow lens in the next few steps.

5. With the hollow lens in place, and the box filled with air, adjust the indices of refraction of each portion. For example, we just did 1: Air, 2: Air, 3: Air. Use the pipette to fill any combination of sections with water or air and make predictions as to whether the light will be refracted or not and whether it will be convergent or divergent.
6. Once you have completed all combinations, fill the box with water and run the same trials. Does anything change?

Lens Surrounded by	Section 1	Section 2	Section 3	Prediction (Converging or Diverging)	Observation (Converging or Diverging)
Air	Air	Air	Air	Neither	Neither
	Air	Air	Water	Converging	Converging
	Air	Water	Air	Converging	Converging
	Water	Air	Air	Diverging	Diverging
	Air	Water	Water	Converging	Converging
	Water	Air	Water	Neither	Diverging
	Water	Water	Air	Neither	Neither
	Water	Water	Water	Converging	Converging
Water	Air	Air	Air	Diverging	Diverging
	Air	Air	Water	Neither	Neither
	Air	Water	Air	Converge	Converge
	Water	Air	Air	Diverge	Diverge
	Air	Water	Water	Converge	Converge
	Water	Air	Water	Convergent	Convergent
	Water	Water	Air	Divergent	Divergent
	Water	Water	Water	Neither	Neither

7. Take one trial into consideration when submerged in water. 1: Water, 2: Air, 3: Air. In this case, what is the “lens” made from
8. Can a Convex Lens be divergent? Can a Concave Lens be convergent? If so, how?

The “lens” is made from air. A convex lens can be divergent when the light passes through the surrounding faster than it passes through the lens. A convex lens can be convergent when the light passes through the surrounding slower than it passes through the lens.

9. Before you leave, please make sure that all the water is put back into your cups and that the lab station is cleaned with a paper towel if necessary.